

**INSTRUCTIONS FOR USE**  
**ZILBRYSQ® (ZIL-brisk)**  
**(zilucoplan)**  
**injection, for subcutaneous use**  
**Single-Dose Prefilled Syringe**

**This Instructions for Use contains information on how to inject ZILBRYSQ.**

**Understanding Your ZILBRYSQ Prefilled Syringe**

Read this Instructions for Use before using the ZILBRYSQ prefilled syringe. This information does not take the place of talking to your healthcare provider about your medical condition or treatment. See the Medication Guide that comes with ZILBRYSQ for important information.

Your healthcare provider should show you how to prepare and inject ZILBRYSQ prefilled syringe properly before you use it for the first time. Do not inject yourself or someone else until you have been shown how to inject ZILBRYSQ correctly.

Keep this Instructions for Use and refer to it as needed until you have used all of the ZILBRYSQ prefilled syringes in the packaging.

For general questions or help, please call your healthcare provider or UCBCares at 1-844-599-CARE (2273).


**How Should I Store ZILBRYSQ Prefilled Syringes?**

- Store ZILBRYSQ prefilled syringes in a refrigerator between 36°F to 46°F (2°C to 8°C) until the expiration date on the carton. **Do not freeze ZILBRYSQ.**
- ZILBRYSQ prefilled syringes may be stored at room temperature up to 86°F (30°C) in the original carton for a single period of up to 3 months. Write the date you removed the prefilled syringes from the refrigerator in the space provided on the carton. After ZILBRYSQ prefilled syringes have been stored at room temperature, do not place them back in the refrigerator. Throw them away if not used within 3 months or if the expiration date has passed, whichever occurs first.
- Keep the ZILBRYSQ prefilled syringes in the original carton before use.
- **Keep ZILBRYSQ prefilled syringes and all medicines out of the reach of children.**

**Important Information You Need to Know Before You Inject ZILBRYSQ**

- **Do not** use ZILBRYSQ if the expiration date on the packaging has passed or the carton seals have been broken.
- **Do not** reuse the ZILBRYSQ prefilled syringe. The prefilled syringe is for 1-time (single use) only. You may get an infection.
- **Do not** inject ZILBRYSQ more than 1 time per day.
- **Do not** miss any doses of ZILBRYSQ. If you miss your ZILBRYSQ dose, inject a dose as soon as possible. Then, inject the next dose at your scheduled time. Do not inject more than 1 dose each day.
- **Do not** use the ZILBRYSQ prefilled syringe if it has been dropped.
- **Do not** remove the needle cap from the ZILBRYSQ prefilled syringe until you are ready to inject.
- **Do not** insert the needle into the skin more than 1 time because this may bend or break the needle, causing trauma to the tissue.
- **Do not** pull back on the ZILBRYSQ prefilled syringe plunger head at any time because this can break the prefilled syringe.
- **Do not** touch the needle guard activation clips at any time because this can cause the early activation of the needle guard.

## ZILBRYSQ Prefilled Syringe Guide to Parts (Figure A):

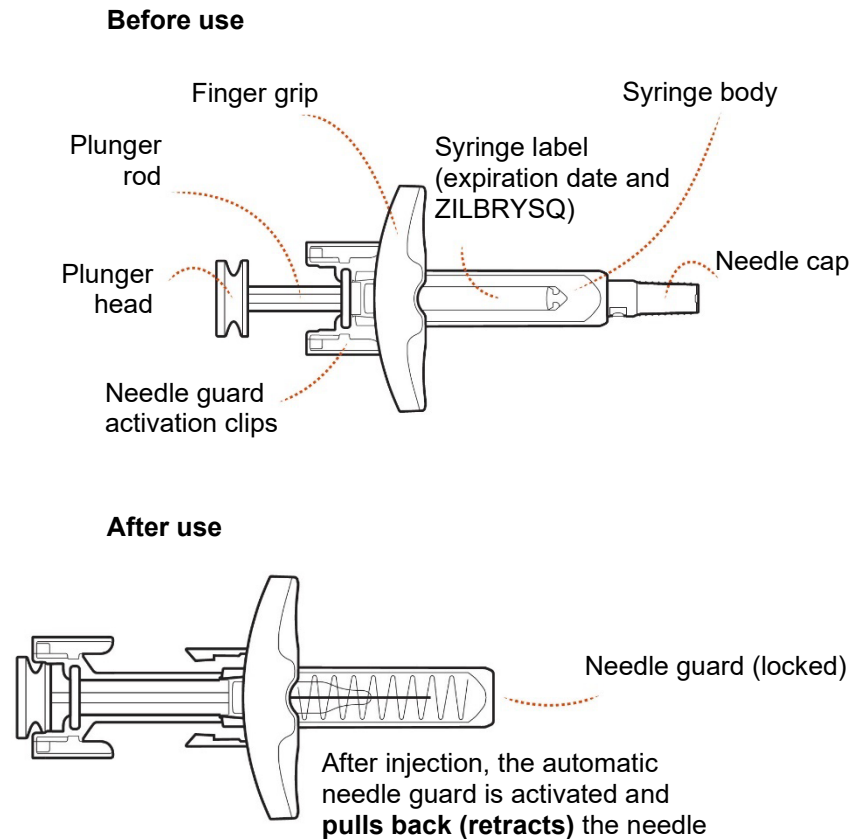

**Figure A**

## Preparing for Injection


### Step 1: Take out the ZILBRYSQ prefilled syringe

#### If the ZILBRYSQ prefilled syringes are stored in the refrigerator:


- Take the ZILBRYSQ prefilled syringes carton out of the refrigerator and remove 1 prefilled syringe from the carton. **Grasp the prefilled syringe body.** Carefully lift the prefilled syringe straight up out of the tray (**Figures B.1 and B.2**). Put the rest of the prefilled syringes in the carton back in the refrigerator.
- Before you inject ZILBRYSQ, let the prefilled syringe warm up to room temperature on a clean flat surface **for 30 to 45 minutes**. This will help to reduce discomfort when injecting.
- Do not warm the ZILBRYSQ prefilled syringe in any other way (for example in a microwave, in hot water, or in direct sunlight).
- Do not remove the needle cap from the prefilled syringe until you are ready to inject.

#### If the ZILBRYSQ prefilled syringes are stored at room temperature:

- Remove 1 ZILBRYSQ prefilled syringe from the carton. **Grasp the prefilled syringe body.** Carefully lift the prefilled syringe straight up out of the tray (**Figures B.1 and B.2**). Any remaining prefilled syringes in the carton should not be placed in the refrigerator after it has been stored at room temperature.

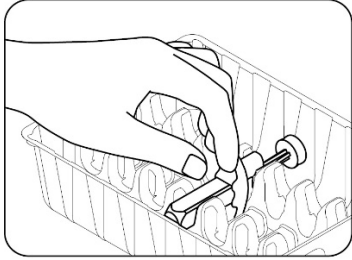


**Figure B.1**

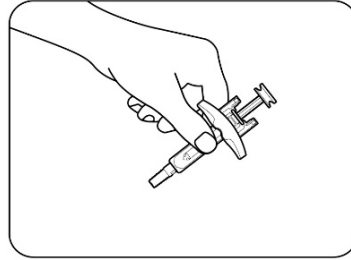

**Figure B.2**

### Step 2: Inspect the ZILBRYSQ prefilled syringe

- Check the prefilled syringe for damage (**Figure C**).
- **Do not** use if any part of the prefilled syringe appears to be cracked, leaking, or broken.
- **Check that** the needle cap is not cracked or broken and attached to the prefilled syringe (**Figure C**).
- **Do not** use if the needle cap is missing or not securely attached.
- Check the expiration date and medicine name (ZILBRYSQ) on the prefilled syringe label (**Figure C**).
- **Do not** use if the expiration date printed on the prefilled syringe has passed.
- **Do not** use if the word ZILBRYSQ does not appear on the prefilled syringe.
- Check the dose appearing on the label (**Figure C**) and make sure it matches your prescribed dose.
- **Do not** use if the dose does not match your prescribed dose.
- Check the medicine inside the prefilled syringe (**Figure C**). The medicine should be clear to almost clear and colorless. It is normal to see air bubbles in the prefilled syringe.
- **Do not** use if the medicine is cloudy, discolored, or contains floating particles.

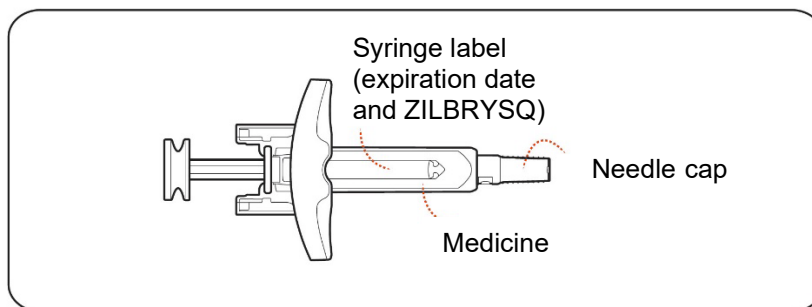


**Figure C**

### Step 3: Gather supplies

Wash your hands with soap and water and dry them with a clean towel.

Gather the following supplies on a clean, flat surface (**Figure D**):

- 1 ZILBRYSQ prefilled syringe
- 1 alcohol wipe (not supplied)
- 1 cotton ball or gauze pad (not supplied)
- 1 adhesive bandage (not supplied)
- 1 sharps disposal or puncture-resistant container (not supplied). See Step 12 for instructions on throwing away the syringe.

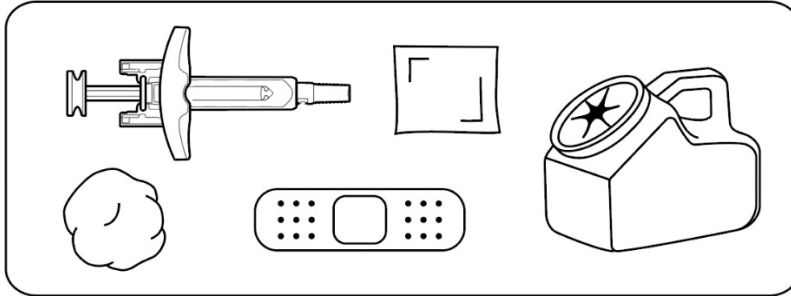


**Figure D**

#### **Step 4: Choose your injection site**


Choose an injection site from the following areas (**Figures E.1 and E.2**):

- The stomach (abdomen), except for the 2-inch area around the belly button (navel) (**Figure E.1**)
- The front of the thighs (**Figure E.1**)

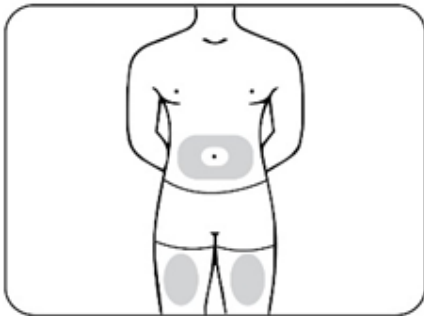

**Figure E.1 – abdomen and thighs**

- The back of the upper arms (only if someone else is giving you the injection) (**Figure E.2**)

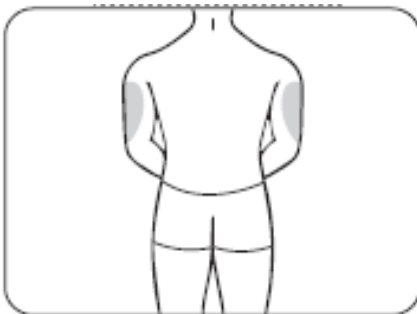

**Figure E.2 – upper arms**

Choose a different site each time you give yourself an injection. If you want to use the same injection site, make sure it is at least 1-inch from a spot you used before.


**Do not** inject ZILBRYSQ into an area that is red, tender, bruised, swollen, hard or that has scars or stretch marks.

#### **Step 5: Clean your injection site**

Clean the injection site using an alcohol wipe (**Figure F**).

Let the skin dry for 10 seconds before injecting.

**Do not** touch the injection site again before giving your injection.

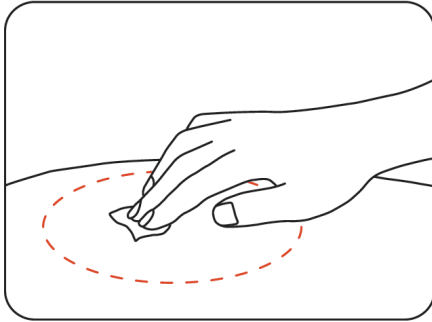


**Figure F**

### **Injecting the Medicine**

#### **Step 6: Remove the needle cap**

Hold the body of the ZILBRYSQ prefilled syringe with one hand and pull the needle cap straight off with your other hand (**Figure G**).

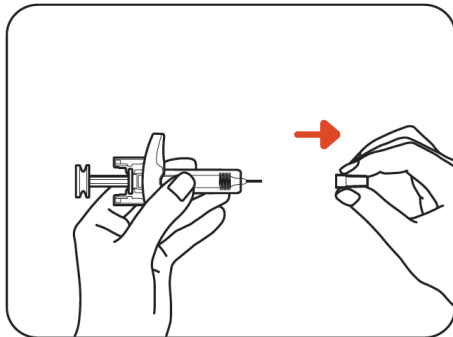

**Figure G**

**Throw away (discard) the needle cap into your household trash or a sharps disposal container.**


**Do not** touch the needle or let it touch anything.

**Do not** recap the needle at any time to avoid injury.

**Do not** try to remove any air bubbles from the ZILBRYSQ prefilled syringe. Air bubbles will not affect your dose and will not harm you. This is normal. You can continue to take your injection.

#### **Step 7: Pinch your injection site**

Use your other hand to pinch the area of cleaned skin and hold it firmly (**Figure H**).

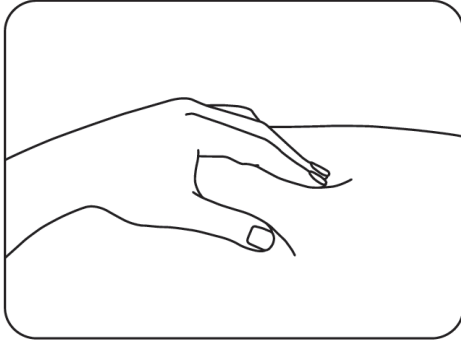


**Figure H**

**Step 8: Insert the needle**

Insert the entire needle into the pinched skin at a 45° to 90° angle. When the needle is fully inserted, hold the ZILBRYSQ prefilled syringe in place (**Figure I**).

**Do not** pull back on the plunger head at any time because this could break the prefilled syringe.


**Do not** touch the needle guard activation clips.

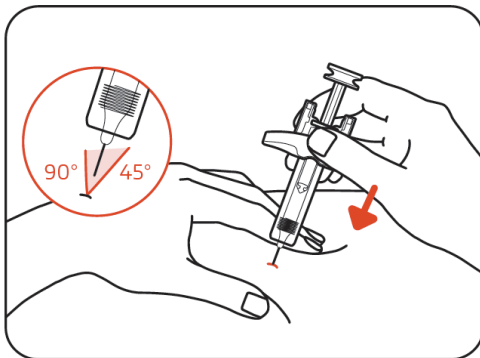

**Figure I**

**Step 9: Release the pinched skin**

When the needle is fully inserted, hold the ZILBRYSQ prefilled syringe in place and release the pinched skin (**Figure J**).

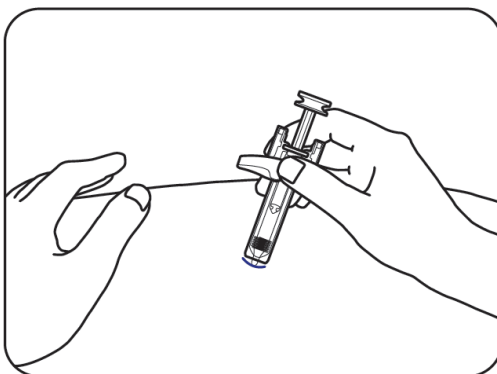


**Figure J**

**Do not** reinsert the needle into the skin if the needle is pulled out when releasing the skin. If this happens, safely throw away (dispose of) the syringe in a sharps disposal container and get a new ZILBRYSQ prefilled syringe to give the injection.

#### **Step 10: Inject the medicine**


Push the plunger head all the way down while holding onto the finger grip to inject all the medicine (**Figure K**).

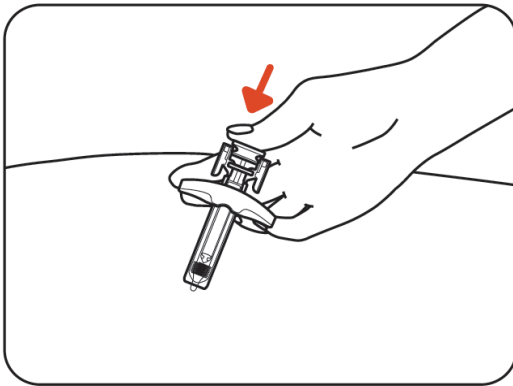

**Figure K**

#### **Step 11: Release the plunger head**

Slowly release the plunger head by lifting your thumb. After a complete injection, the needle guard will cover the needle and you may hear a click (**Figure L**).

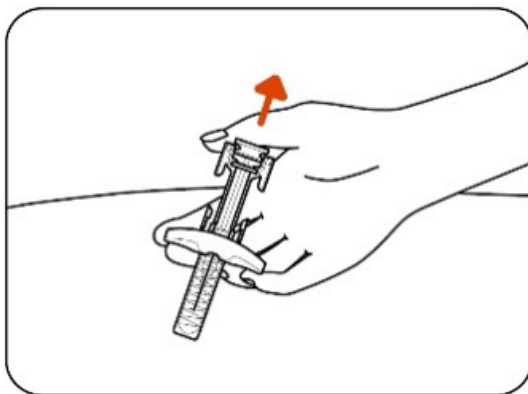

**Figure L**

#### **Step 12: Throw away (dispose of) the used ZILBRYSQ prefilled syringe**

Throw away (dispose of) the used ZILBRYSQ prefilled syringe into a sharps disposal container (**Figure M**) right away.

**Do not** throw away (dispose of) the ZILBRYSQ prefilled syringe in your household trash.

If you do not have an FDA-cleared sharps disposal container, you may use a household container that is:


- made of a heavy-duty plastic,
- can be closed with a tight-fitting, puncture-resistant lid, without sharps being able to come out,
- upright and stable during use,
- leak-resistant, and

- properly labelled to warn of hazardous waste inside the container.

**Do not** throw away (dispose of) your used sharps disposal container in your household trash unless your community guidelines permit this.

**Do not** recycle your used sharps disposal container.

Always keep the sharps disposal container out of the reach of children.

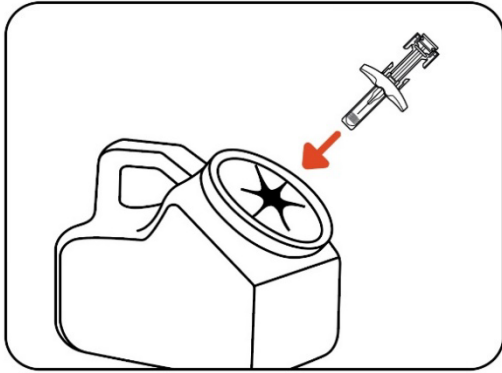


**Figure M**

**Note:** When your sharps disposal container is almost full, you will need to follow your community guidelines for the right way to dispose of your sharps disposal container. There may be state or local laws about how you should throw away your used needles and syringes. For more information about safe sharps disposal, and for specific information about sharps disposal in the state you live in, go to the FDA's website at: <http://www.fda.gov/safesharpsdisposal>.


**Step 13: Examine the injection site**

Press a cotton ball or gauze pad over the injection site and hold it for 10 seconds (**Figure N**).

**Do not** rub the injection site. You may have slight bleeding. This is normal. Apply an adhesive bandage if needed (**Figure O**).

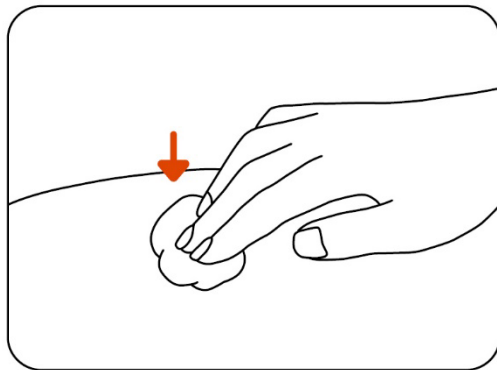

**Figure N**

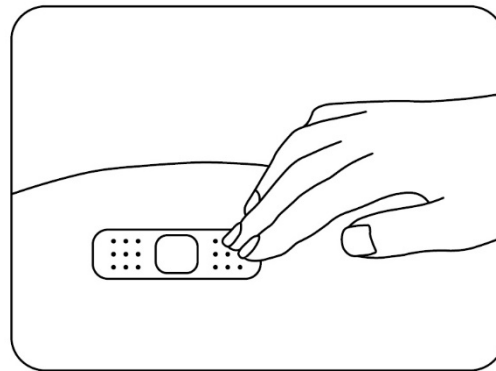

**Figure O**

**The injection is complete.**

For more information or help, contact your healthcare provider.

**If you have any questions, call UCBCares at 1-844-599-CARE (2273) or visit [www.AskUCBCares.com](http://www.AskUCBCares.com).**

Manufactured for:  
UCB, Inc.  
Smyrna, GA 30080

This Instructions for Use has been approved by the U.S. Food and Drug Administration.

Issued: 10/2023